spirax sarco

**TI-P180-11** CMGT Issue 21

# BTM7 and BTM7 Food+ Balanced Pressure Thermostatic Clean Steam Trap


# **Description**

#### BTM7

The Spirax Sarco BTM7 is a maintainable thermostatic steam trap designed to remove condensate from clean (chemical free) steam systems with minimal condensate retention. Applications include steam barriers, block and bleed installations, mains drainage and CIP/SIP of vessels and reactors and process lines. Manufactured in 316L stainless steel with minimal crevises and a typical internal surface finish of 1.6 - 3.2 µm Ra, it is self-draining and operates close to steam temperature. Each trap is individually packaged within an "ISO CLASS 7" clean environment with protective end caps and sealed in a protective plastic bag.

#### BTM7 Food+

Designed, manufactured and approved for Steam and Condensate applications. The BTM7 Food+ range complies with EC1935:2004 Food Contact Materials. It also complies with regulation EC2023:2006 on good manufacturing practice for materials and articles intended to come into contact with food. A list of the materials that could come directly or indirectly into contact with foodstuffs can be found in the Declaration of Compliance available for this product.

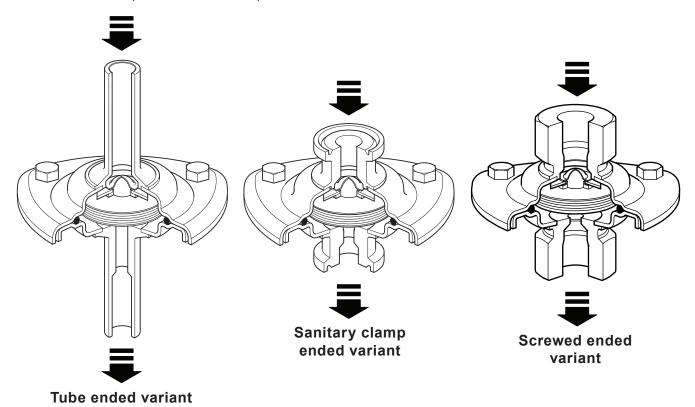

## Options and available types

Contact Spirax Sarco for further information

Fixed bleed to ensure fail open operation.

Special connections to suit most piping systems.

## **Standards**

- The BTM7 has been designed and manufactured in general accordance with ASME BPE standards.
- The unit also complies with the requirements of the European Pressure Equipment Directive 2014/68/EU.
- FDA CFR Title 21. Paragraph 177.2600.
- USP Class VI Biological Reactivity Cytotoxicity Testing In-Vivo <88> extracted at 121 °C for 1 hour.
   ADI Free (Animal Derived Ingredients) for materials used, manufacturing processes involved in producing the part.
- TSE/BSE Free as Certified.
- Full material lot number product traceability.

## Sizes and pipe connections

	Standard	Issue	Туре	Sizes						
	A CME DDE		Type A			1/2"	3/4"			
	ASME BPE		Type B					1"		
Sanitary clamp ended variants	DIN32676	2001-02	Series 2			DN15				
		2009-05	Series A			DN15				
			Series B	DN8	DN10	DN15	DN20	DN25	DN32	DN40
			Series C			1/2"	3/4"	1"		
	DIN11866	2016-11	Series A			DN15				
			Series C			DN15	DN20	DN25		
Tube ended variants	DINI44050	1999-01	Series 1		DN10	DN15				
	DIN11850		Series 2			DN15				
	ISO1127	1997	Series 1	DN8	DN10	DN15				
Screw ended variants	Socket BSP T Rp (ISO 7-1) Socket NPT				1/4"	1/2"	3/4"	1"		

#### Notes:

- 1. On request other connection options are available at extra cost. Please note that seat end spares for specially requested connections will require a minimum order quantity Please consult Spirax Sarco for further information.
- 2. BTM7 Food+ available as tabulated above for ASME BPE and BSP connections. NPT available in 1/4" ONLY.


## Certification


This product is available with the following certification:

- EN 10204 3.1 Material Certification Pressure Retaining Parts Chargeable
- EN 10204 3.1 Material Certification for Wetted Parts including the WFI element fill (also available for element spares) Chargeable.
- Typical Internal Surface Finish F.O.C.
- Certificate of Compliance for FDA, USP Class VI Testing Statement, and ADI Free Statement F.O.C.
- TSE-BSE Statement F.O.C.
- EC1935:2004 Declaration of Compliance F.O.C. (Food+ range only)
- (EC)2023:2006 Good Manufacturing Practice for Materials and Articles Intended to come into Contact with Food.
- Declaration of Conformity BS EN ISO 14644-1:2015 Class 7 Clean Room F.O.C.
- Typical Test Report F.O.C.


**Note:** All certification/inspection requirements must be stated at the time of order placement, and may be subject to additional cost as detailed above.

# **Materials**

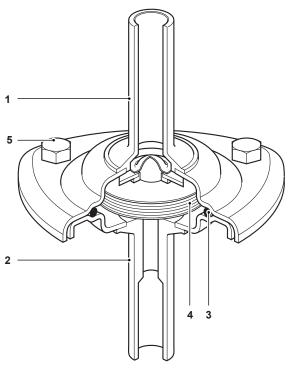

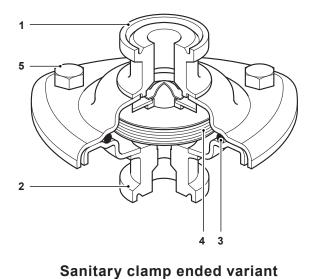


**Tube ended variant** 

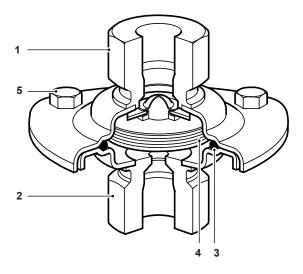

Screwed ended variant

No.	Part	Material	
1	Body (inlet)	Stainless steel	316L (1.4404)
2	Body with seal (outlet)	Stainless steel	316L (1.4404)
3	'O' ring	FKM	
4	Capsule element	Stainless steel	316L (1.4404)
_	Nuts and bolts	Stainless steel	BS 6105 Gr. A4 80
5	Washers	Austenitic stainless steel	


# Pressure/temperature limits

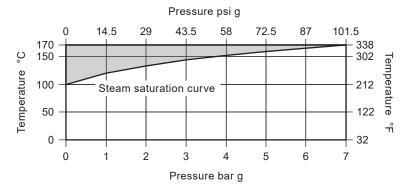


The product should not be used in this region as damage to the internals may occur.


Body	design conditions		PN7	
PMA	Maximum allowable pressure	7 bar g @ 170 °C	101.5 psi g @ 338 °F	
TMA	Maximum allowable temperature	170 °C @ 7 bar g	338 °F @	
Minimum allowable temperature		-10 °C	14 °I	
РМО	Maximum operating pressure for saturated steam service	7 bar g	101.5 psi g	
ТМО	Maximum operating temperature	170 °C	338 °F	
Minimum operating temperature		0 °C	32 °F	
Designed for a maximum cold hydraulic test pressure of:		10.7 bar g	155.1 psi g	


# **Capacities**

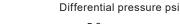

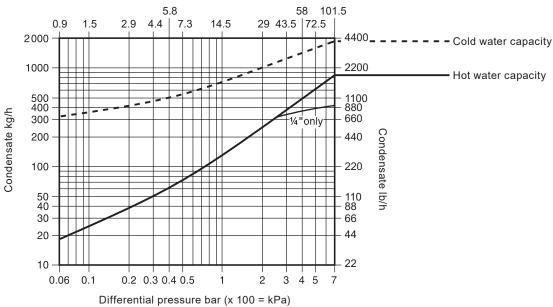


# Dimensions/weights (approximate) in mm (inches) and kg (lbs)

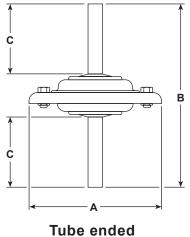

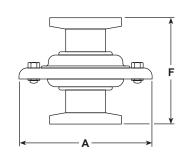

variant

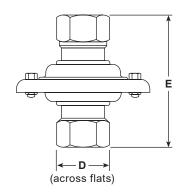

Sanitary clamp ended variant

Screwed ended variant

Products each weigh 1 kg (2.2 lbs) ma							maximum '	
Standard	Sizes	A	B Tube	С	D Screwed	E Screwed	F Sanitary Clamp	Weight
Sanitary Clamp ended va	riants							
ASME BPE Types A and B, DIN32676 Series 2 + A, B, and C Variants	See Sizes and pipe connections table	71 (2.8)					49 (1.9) for DN8 to DN20, and 49 (1.9) for	1 kg * (2.2)
							DN25 to DN40	
Tube ended variants								
DIN11866 Series A and C, DIN11850 Series 1 and 2, ISO1127 Series 1	See Sizes and pipe connections table	71 (2.8)	107 (4.2)	38 (1.5)				1 kg * (2.2)
Screw ended variants								
Socket BSP + NPT	See Sizes and pipe connections table	71 (2.8)			27 (1.0) 27 (1.0) 32 (1.2) 41 (1.6)	58 (2.3) 74 (2.9) 81 (3.2) 95 (3.7)		1 kg * (2.2)

## Safety information, installation and maintenance

For full details see the Installation and Maintenance Instructions (IM-P180-05) supplied with the product.

#### Installation note:

The trap is designed for installation in vertical lines with the flow downward to ensure self-draining operation. Do not expose element to superheat conditions since over-expansion may result.

Suitable isolation valves must be installed to allow for safe maintenance/replacement.

# How to order

**Example:** 1 off Spirax Sarco DN15 BTM7 maintainable thermostatic clean steam trap with tube ends to ISO 1127, Series 1. (21.3 mm O/D x 1.6 mm wall thickness). Tangent length of tube ends to be 40 mm for ease of orbital welding. Body to be self-draining. Suitable for pressures up to 7 bar g.

## Spare parts

Available parts are shown in solid outline. Parts drawn in a grey line are not available as spares.


#### Available spares

Capsule element assembly	4
'O' ring (packet of 3)	3
Body with seat (outlet) - state connections	2

#### How to order spares

Always order spares by using the description given in the column headed 'Available spares' and state the size, type and end connection of the trap.

Example: 1 off Element assembly for a Spirax Sarco 1/2" BTM7 thermostatic clean steam trap having screwed NPT connections.

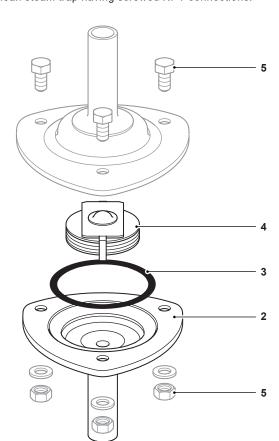

# Recommended tightening torques

Item	Part	or mm	<b>*</b>	N m	(ft lbs)
5	Nuts and bolts	8 A/F	M5	3-4	(2.2-3)